




FMT CLASSIC

FLOOR-MOUNTED TUBE STAND SYSTEM

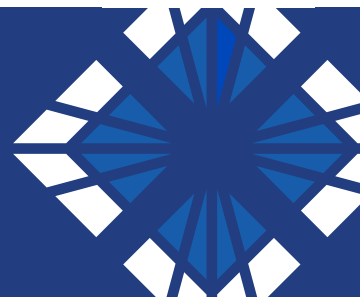

Maximum Capability, Low Cost of Ownership

The **FMT Classic** X-ray system delivers maximum imaging capability with an economical and compact design. The low ceiling height requirement of 92.5" maximizes your existing space, avoiding costly room modifications.

Plus the FMT Classic comes with the industries best **5-year parts warranty** making it the ideal solution for medium to high-volume clinics that have limited space and are looking for a cost effective system.



Technologist and Patient Friendly Features

The tube stand allows for patient imaging on the table or wall stand and can also accommodate wheelchair exams. Technologist-friendly features, such as the easy-read angulation dial and one-touch release handle make positioning easy and helps to improve patient throughput.

Safety features like the recessed foot pedals and positioning features like the comfortable patient chin rest and table working height indicator allow the technologist to focus on the patient, improving the quality of the patient experience.

The FMT Classic system includes:

J700 FLOOR-MOUNTED TUBE STAND

Transverse and rotational movements with freestanding installation capability.

- ✓ 40" and 72" electronic detents for quick horizontal positioning for all weight bearing and upright exams
- ✓ The easy-read backlit blue angulation dial makes identifying and adjusting angles fast and accurate



S211 RADIOGRAPHIC TABLE

The elevating, four-way float, smooth top allows for easy transfer and positioning of patients.

- ✓ 400 pound patient capacity
- ✓ Recessed four pedal design provides easy operation and patient safety
- ✓ Indicator light for optimal working table height

S109 WALL STAND

Wall stand with rigid construction and electric locks ensures precise alignment.

- ✓ Fully counterbalanced center-mount
- ✓ Push-button lock release handle for easy positioning
- ✓ Comfortable patient chin rest

FMT Classic Specifications

Radiographic Table (S211)

ELEVATING FOUR-WAY FLOAT TOP TABLE

- 84 in. x 30 in. flat, seamless table top
- Variable table height, 23 in. to 35 in.
- +/- 5 in. transverse travel
- 31 in. longitudinal travel
- Foot treadle operation of table top and height movements
- Grid cabinet with electric lock
- Grid: 10:1 ratio, 103 LPI
- Heavy-duty deluxe cassette tray 06338
- Operates on 120VAC
- 400 lb. rated patient load

Wall Stand (S109)

WALL-MOUNTED WALL STAND (17"X 17")

- Fully counterbalanced
- 17 in. x 17 in. grid cabinet with front panel travels from floor to 76 in.
- Grid: 10:1 ratio, 103 LPI, specify focal range
- Specify left or right load
- Heavy-duty deluxe cassette tray 06338
- Electric locks
- Operates on 24VDC

TUBE STAND (J700)

FLOOR-MOUNTED TUBE STAND

- 10 ft. rails
- Rotational base, 360° with 90° detents
- Transverse tube arm with 8.5 in. travel
- Angulation dial, operator handgrips and electric locks
- Platform tube mount

Generators

- Choose from our 40 kW - 80 kW HF generator options



SUMMIT INDUSTRIES
A NETWORK OF SERVICE & SUPPORT

Summit Industries' successful service network surrounds you with a dynamic team of exceptional experts. Your organization and your patients will benefit from the experience and best practices that only Summit Industries and their Authorized Dealer Network can offer.



SUMMIT INDUSTRIES, LLC

7555 N. Caldwell Avenue, Niles, Illinois 60714
(p) 773-353-4000 | (p) 1-800-729-9729



LITR031-0225



www.amradmedical.com